

CALMERA DELUXE – Massage Chair



Product user manual

SL Super Long SL Track One
Button Sliding Forward



Living Art

Thank you for your purchase of our massager.

This massage chair adopts the ergonomics design to realize the multi parts massage on the back, the waist, the hip, the leg, the arm and so on. It is a luxury massage chair with many kinds of massage, shiatsu, kneading, flapping, knocking,rolling, air squeezing and so on. Please read this manual carefully before using it, and pay special attention to safety precautions to ensure proper operation and use of this product. Keep the user's handbook in good order for future enquiry.



CATALOG

Security matters	03
Product structure and parameters	06
Remote control operation instructions	07
Instructions for the operation of the massage procedure	08
Installation steps / methods	09
Daily maintenance and troubleshooting	10



Security matters

Please read all the instructions in detail before using this product.

In order to reduce the risk of fire, burns, electric shock or injury to others, follow the following instructions:

- 1) When no machine is used, please make sure that the electrical plug of the electrical server is removed from the socket. It is necessary to remove the plug from the socket in a timely manner without the use of the machine or the device / disassembly parts.
- 2) To avoid fire, electric shock or injury caused by overheating, do not place products under other covers such as pillows or blankets.
- 3) If the product is used or close to children, physical and mental handicapped and disabled, adult guardians should be supervised and directed.
- 4) Be sure to use this product in accordance with the purpose described in the manual. Do not use the spare parts recommended by the non manufacturer.
- 5) When the following condition occurs, please do not use this product: the power cord or the power hand plug is damaged, the operation is abnormal, the ground or other damage has been dropped, and the water has been dropped. Please send the product back to the maintenance service center for inspection and maintenance.
- 6) Do not drag the power cord to move the product, or use the power cord as a handle.
- 7) The power cord must be kept at a safe distance from the heated surface.
- 8) Do not operate the product under the condition of air outlet obstruction. Please keep the outlet smooth and avoid clogging with cotton dust and hair.
- 9) Do not insert or drop objects into any openings of the product.
- 10) Please do not use this product outside.
- 11) Please do not operate the product in the place of spray product or oxygen filled place.
- 12) Cut off the power supply. Please close all the operation control and remove the power plug from the socket.
- 13) Be sure to supervise and confirm that children do not use massage chairs as amusement equipment to ensure their safety.



1、 Security matters

- Do not use massages in excessive humidity or dust dense environment, so as to avoid electric shock or machine failure.
- Do not use massages at room temperature above 40°C.
- Do not use the massager as near the electric heater, near the stove or under the direct sunlight.
- The massager is designed for personal use and is not suitable for any commercial purpose.

Be careful

- The massager is limited to indoor use.
- If the massager is used in low chamber temperature, do not suddenly raise the temperature.
It is recommended that the temperature be gradually adjusted to normal range before use.
- If the massager is moved from a cryogenic storage to a warm environment, it may affect the normal operation of the machine and even the failure of the machine because of the accumulation of water vapor in the machine parts.

2、 Safe use

- To ensure safe and proper use of massages, please read the operation manual carefully before using them.
- Do not use if the sheath of the massager has been damaged.
- To avoid reducing or injure the massage effect, do not use the massager or other electric appliances such as electric blanket.
- Do not use a massager within one hour after eating, so as not to cause any adverse effect. Do not use this product when your hands or other parts of the body are moist.

- Please do not insert anything or flammable goods into the massager. This can lead to electric shock, malfunction or fire.
- When using, do not cover the massager, so as not to cause too much trouble or fire.
- If your body feels any discomfort during the massage, please stop immediately.
- This product is not used for self diagnosis or treatment: if it is needed for diagnosis or treatment, it should be handled by professional medical personnel.
- If you are receiving medication or medical conditions, consult your doctor before using it. If you feel unwell or are pregnant, please do not use this product.
- When you are pregnant or in the physiological period, consult the doctor first and follow the doctor's advice.
- It is suggested that the time limit is 15~30 minutes. Before using it again, please confirm that the massager has been unpowered and stopped for more than 10 minutes. If continuous use for a long time, it may start overheating protection device, automatic power failure.
- People who are not sensitive to heat on surfaces of appliances must be noted when using them.

Be careful

- Please avoid using a massager after drinking or in untimely use.
- Please do not directly massage the head, elbow, knee, chest, abdomen and so on.
- When inadvertently sprinkle or dump the liquid on the massager, please remove the power plug from the socket immediately. When the hands are wet, do not plug or remove the power plug.
- Please do not damage the power cord or operate the massager under the condition that the power cord is damaged. This may cause an accident of electric shock or machine failure.
- If the following conditions occur, please turn off the switch immediately, cut off the power supply, and contact our company.
 - Accidentally splashing water on the massager may cause electric shock.
 - The cloth sleeve of the massager is damaged and the internal parts are exposed.
 - If there is abnormal or severe pain during massage, please stop the massage immediately and seek medical advice.
 - If you find any malfunction or abnormal condition during use, please stop using it immediately.
 - In case of interruption in use, care should be taken to avoid injury to users due to sudden recovery of power.
 - When thunder and lightning occur.

Product structure and parameters




Product Name	Deluxe Massage Chair
Model	
Product Power	135 W
Input Voltage	220V—240V
Frequency	50~60Hz
Working Hours	20 min

Remote control operation instructions



Shortcut

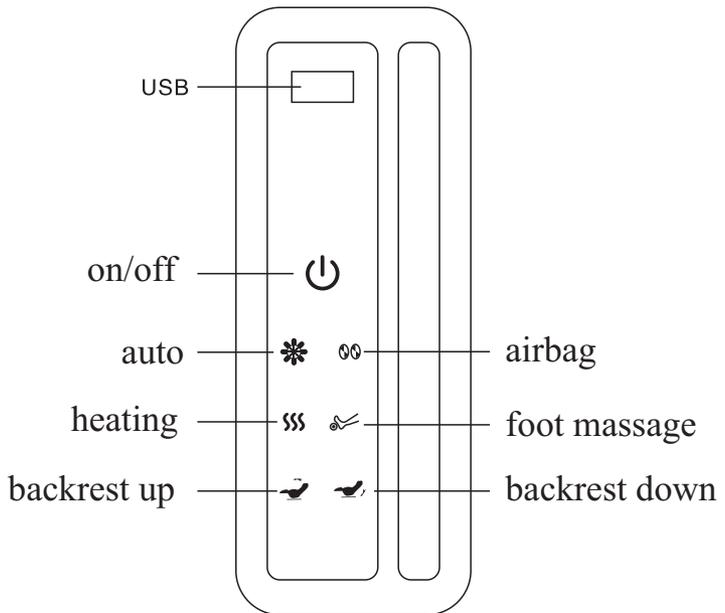

Instructions for the operation of the massage procedure

1. Standby

[Standby]: Press the “” button and the whole machine enters the automatic massage working state, and then press the “” massage chair to stop working.

2. Timing

[Timing]: Press the “” button in the power-on state to select different massage time.

3. Automatic

 [Auto]: By pressing the automatic button: “Auto” select a variety of different automatic massage programs to meet your different needs.

 **Soothing mode**—a pleasant and relaxing full-body massage to relax your mind and body.

 **Powerful message**—relaxing massage to adjust the pressure after work.

 **Stretching mode**—relax muscles and vitalize bones, restore muscles and bones vitality.

 **Sleep mode**—improve sleep and promote physical and mental health.

 **Shoulder mode**—relax the shoulder and neck muscles.

 **Waist mode**—relax the back muscles.

 **Ladies mode**—relax the full body

 **Elderly mode**—rejuvenate the full body

4. Manual

[Manual]: Press the “M” button in the power-on state. Turn on specific massage functions such as kneading, beating, kneading, tapping, and acupressure.

5. Fixed point

[Fixed point]: In manual mode, press the “” key to start the fixed point positioning massage.

6. Shoulder height

[Shoulder Height]: Press the “” button under the power-on state, the backrest motor returns to the shoulder position, and press the up and down buttons to adjust the best position of the movement.

7. Airbag automatic

[Airbag Auto]: Press the “” button in the power-on state to open or close the whole body airbag.

8. Air pressure adjustment

[Air pressure adjustment]: Press the “” button in the power-on state to turn on the function of air bag massage intensity adjustment.

9. Massage speed

[Massage speed]: Press the “” button in the power-on state to adjust the massage intensity.

10. Heating

[Heating]: Press the “” button in the power-on state to turn on or off the heating function.

11. Pause

[Pause]: Press the “” button in the power-on state to stop the machine.

12. Width adjust

[Wide Tuning]: Press the “” button in manual mode to adjust the width, middle and narrow of the massage head spacing

13. Zero gravity

[Zero Gravity]: Press the “” button in the boot state to turn on the zero gravity massage mode.

14. Roller

[Roller]: Press the “” button in the boot state to turn on and off the roller massage function.

15. Backrest up /backrest down

[Backrest up]: Press the “” button in the power-on state, and the chair back will rise up.


[Backrest down]: Press the “” button in the power-on state to lie on the back of the chair.

16. Calf up / Calf down

[Leg Raise]: Press the “” button in the boot state to control the raising of the calf.

[Leg Down]: Press the “” button in the power-on state to control the lowering of the calf.

Installation steps / methods

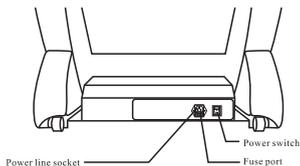

When you move, press the back and raise the front

Attention

- ① Please check and confirm if the power cord has been pressed under the massage chair.
- ② Whether the ground line is connected well
- ③ Before inserting the plug, make sure that all the operation buttons, including the power switch, are turned off.

Daily maintenance and troubleshooting

- 1, Daily cleaning, please use neutral detergent and soft cloth wipe, do not use corrosive detergent.
- 2, If the power cord is damaged, there must be a manufacturer's replacement for the maintenance shop to avoid danger.
- 3, When you are not using this product or cleaning, please remove the power source to avoid danger.
- 4, Do not use any sharp items to cut the surface of the product.
- 5, Do not put products in damp environment, should be placed in dry and ventilated place.

Failure phenomenon	Possible cause	Processing method
Sound (running sound) The inflatable sound of the air bag The bleed sound of the air bag Switching sound of solenoid valve Working sound of air pump	When the product works The sound of a mechanical operation	The sound is normal. There is no need to deal with it.
The whole machine is completely out of work	The power plug is not inserted into the socket The socket switch is unopened The massage button is not pressed	Check the power plug and socket, ensure that the plug is in good contact with the socket, press the massage button.
The massage was stopped in the course of a sudden stop	This product has just completed 15 minutes of rated working time and automatically switched to standby mode.	If you want to massage once again, please press the key to open, open massage function
	Overheating protection of motor over rated use time	Shut down, cut off the power supply to motor cooling after use