


SATORI ELITE

Massage Chair

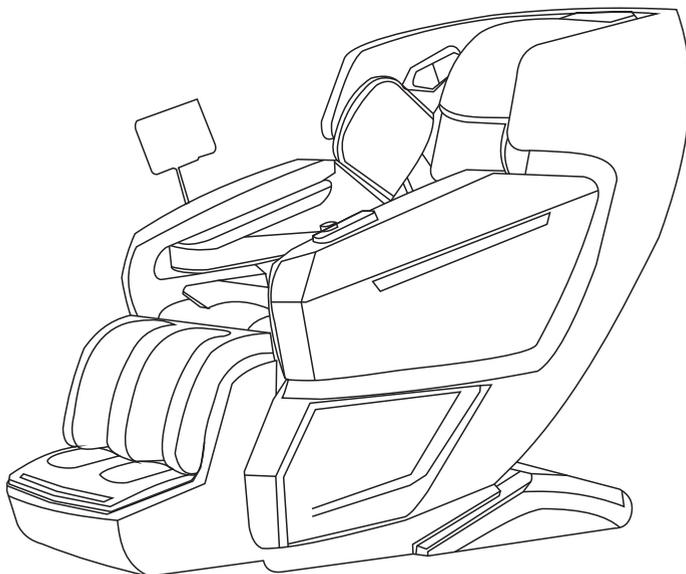

catalogue

ontents

Safety precautions	1-4
Name of each part	5-7
(1) Appearance structure	5
(2) Manual panel	6
Usage method	7-9
Massage chair installation	10-12
(1) About the place of placement	10
(2) The method of moving in and out	11
(3) Installation method	12
Maintenance method	12
Common problems and solutions	13
Product specification	14

Safety Precautions



Safety Precautions

- Please read this [Safety Precautions] carefully before use and use this product correctly.
- The precautions here provide tips on how to use the product safely and correctly to prevent possible harm or injury to the product, the user, and others.
- This precaution is distinguished by two prompts: [Warning] and [Caution] based on the size and severity of the hazard or injury that may result from improper operation. Content related safety is very important, please be sure to abide by it.

 Warning:	Misoperation may cause serious injury to the user.
 Note:	Misoperation may result in user injury and property damage. Physical damage refers to damage related to houses, property, livestock, pets, etc.

【 Picture 】

	This sign indicates warning or caution content The picture on the left is generally used after adding a warning or caution.
	This sign  informs of prohibited content, and the picture shows a description of the specific prohibited content.(The picture on the left shows that disassembly is prohibited)
	This sign informs you of mandatory or required content.

- After reading, please be sure to save it in a place where users can see it at any time.



Warning



- If smoke or abnormal noise occurs (abnormal state), please immediately disconnect the power supply and confirm safety, please contact the dealer for repair. It is strictly prohibited to repair it yourself to avoid danger.



- There is water inside the host.
- Debris fell into the host.
- Heavy fall, fall, or case damage.



Notice



Must be obeyed

- During use, if you feel any abnormality in your body, you should stop using it immediately and consult a doctor.
- This product should be placed on a flat ground for use to avoid injury accidents.
- Before using this product, please confirm whether the power switch is working properly to avoid injury accidents.



Safety precautions

Safety Precautions

Warning

	You should not use this device
	<p>1:Osteoporosis patients. 2: Heart problems and electronic medical devices such as pacemaker 3:Pregnant women or menstruating women. 4: A person with a physical or physical Condition. 5:Fever allergies. 6:The doctor instructs the person to be recuperate or to feel ill.</p>
	<p>This product massage program is repeated for 15 minutes. Please avoid continuous massage in the same part. If you need to use it continuously, please wait for more than 30 minutes to avoid excessive stimulation of muscle or nervous system.</p>

Warning

	Do not use the machine under voltage that is not suitable for the local standard.
	Be sure to use the plug that matches the plug.
	Please use this product in this manual.Please do not use this product outside.
	This appliance is not intended to provide physical disabilities,sensory And neurological
	deficiencies or lack of experience and common sense (including children).Unless they are given the supervision and guidance on how to use the instrument for their safety and responsibility;Children must be supervised to ensure that they do not play with instruments.If the power cord is damaged,it must be replaced by the manufacturer,the maintenance department or similar departments in order to avoid danger.

Safety precautions

Safety Precautions



Warning



The massage chair change the angles automatically during the automatic procedure, and it is forbidden to contact any gaps and gaps to avoid injury.

The power cord must not be wrapped around the product to avoid electric shock, short circuit or fire.

Do not allow children to ride on this product to avoid injury or machine malfunction.

When back, seat and the legs adjusting, make sure there are no other people, pets or items under the rear and leg part, and do not put any part of the body, including the hands and feet, during the adjustment process. Including Into the backrest and seat, backrest and armrests, armrests and shoulders, seat and the cover under the armrest to prevent injury.

Do not use the product naked to avoid skin damage.

Be sure to use the product in the manner described in this manual. Do not sit directly on the legs or armrests to avoid injury or accidents.

Do not use the machine if it is damaged, so as to avoid injury or accident.



Attention



If not use or the power is off during use , please turn off the power switch immediately and unplug the power cord to avoid leakage, electric shock or fire.



Do not touch the power outlet, plug or power cord with wet hands to avoid electric shock.



When cleaning, be sure to turn off the power switch and unplug the power cord from the power outlet to avoid fire or electric shock.



After a meal or drinking, you should not use this product immediately to avoid discomfort.

Do not hang heavy objects on the power cord or manual control line. The hand controller should not be impacted or dropped by external force to avoid malfunction, electric shock and fire.

Do not use this product near fire or heaters, and do not use it on floors with floor heating ..Do not smoke when using, to avoid fire



Safety precautions

Safety Precautions



Attention:



- Do not pull or use a knife, hard sharp objects to scratch the fabric to avoid injury or accident
- Do not bend, pull the power cord, hand control line forcibly, to avoid machine failure or fire accident.
- The hand controller cannot be impacted or dropped by external force to avoid malfunction or machine failure.
- Do not use this product when the power cord is damaged or the power socket is loose, so as to avoid electric shock, short circuit or fire accident.
- Please do not put your fingers into any gaps when the leg lifting or down, to avoid injury.
- Do not force external forces on the components in each action to avoid personal injury or machine malfunction.
- Do not use while asleep, unintentional use may cause injury and physical discomfort.
- Do not drag on the wooden floor, carpet or easily scratched ground in the direction of the wheel.
- Gently place after moving to avoid scratching the ground



Attention



Do not cover a blanket when using this product to avoid a fire accident caused by excessive internal temperature of the machine.

Do not stand on the fuselage or on the leg parts, step on or stack heavy objects to avoid personal injury or machine malfunction.

Do not use this product near fire or flammable materials to avoid product deformation, machine failure or fire accidents.

Do not wipe or spray this product with organic solvents, alcohol or insecticides to avoid cracks or electric shocks.

Do not let children or people with intellectual disabilities approach during use to avoid injury.

People who are not allowed to use this product to avoid discomfort or injury


Name of each part





-Function indication

Appearance structure

▲ This machine has automatic detection of body shape function, please relax and lie down when starting!

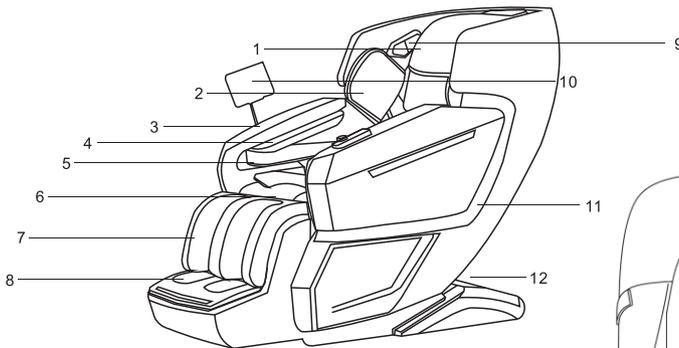

- | | |
|--------------------|-----------------------|
| 1. Headrest | 8. Feet |
| 2. Shoulder clamp | 9. Negative ions |
| 3. Shortcut keys | 10. Manual device |
| 4. Clamp your hand | 11. Atmosphere lights |
| 5. Health module | 12. Base |
| 6. Cushions | |
| 7. Legs | |

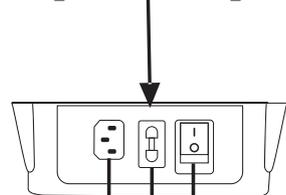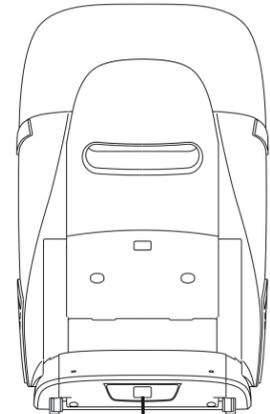

Power line interface

Fuse 3.15A

Power switch




Part name

Function —————

Controller



Controller installation instructions

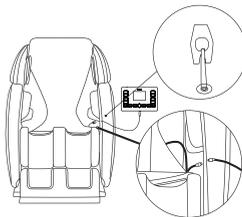


- Place the controller bracket holder on the screw of the left armrest, and tighten the bottom of the hand control bracket with the screw;
- Find the hand control seat cushion connector on the left side, connect the hand control connector and the seat cushion connector, and tighten them;
- Place the hand control on the hand control bracket and the hand control installation is completed.
- Do not pour water or other liquids onto the massage chair or hand control
- Do not operate with wet hands.
- Do not sit on, drop or step on the hand controller, and do not forcibly pull the wire of the hand controller.
- Do not press the liquid display and touch panel with sharp or hard objects. (Otherwise it may cause damage)

Instructions for use



Setting

In order to achieve the best massage experience, an 8-second body shape detection is performed when the device is turned on. For accurate detection, please keep your head against the headrest. The user presses the "up" and "down" keys next to "OK" to perform manual fine-tuning.

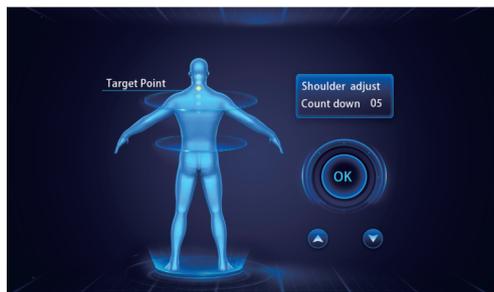

Press the menu key and click Master Smart Selection. Automatic mode: Select the corresponding automatic mode program. (9 modes including sleep mode, neck and shoulder, Thai stretching, lumbar massage, strong massage, elderly mode, queen slimming, relaxing and relaxing, and exercise recovery)



Press the Menu key and click the Home key
Home page content: page settings, power on and off settings, pause settings.





Instructions for use

enitortnarl

Setting

Press the menu key and click the manual key.

Manual mode - Technique: Select the corresponding massage technique.
(Total 20 techniques to choose from)



Press the menu button: click on the up movement button.

Manual mode → Technique: Select the corresponding massage technique.
(20 techniques to choose from, including kneading, knocking, tapping, shiatsu, kneading+knocking, and kneading+tapping, etc)



Press the menu button: click the lower movement button.

Manual mode → Technique: Select the corresponding massage technique.
(20 techniques to choose from, including kneading, knocking, tapping, shiatsu, kneading+knocking, and kneading+tapping, etc)



Instructions for use

Instructions



Setting

Press the menu key: click the sitting posture key

Sitting mode: 1. Chair back lift adjustment, 2. Leg lift adjustment,
3. Zero gravity automatic program.



Press the menu key: click the settings key

Setting mode 1. English voice switching, 2. Bluetooth connection, 3. A i voice



Press the menu key: Click the air bag key.

Air bag mode: You can choose to turn on or off the corresponding air bag.
(Strength can be adjusted)Can choose foot roller,heating turn on or off.






Massage chair installation

Massage chair installation

Placement

Please ensure that there is 35cm in front of the massage chair and about 8cm in the back to avoid hitting the wall or other objects when the massage chair is tilted.

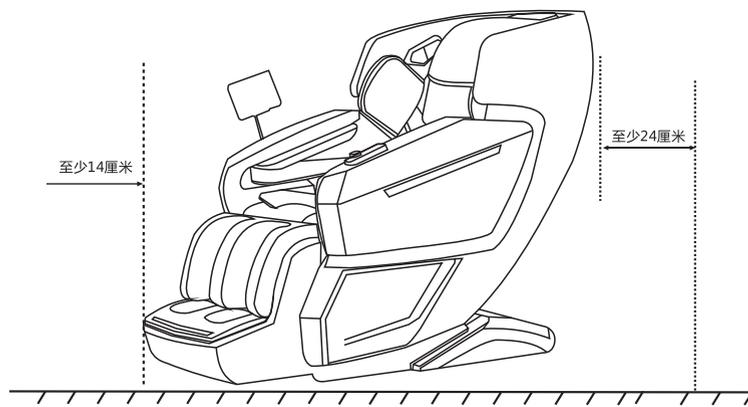

Please leave adequate space around the massage chair when preparing for massage



Notice

- Do not use in places with high humidity such as bathrooms or saunas.
- Please place the machine on a flat and non-slip surface.
- Do not use it in places with a lot of debris, dust, or places with high ambient temperatures (above 40C)
- Do not place the machine in a place exposed to direct sunlight, or use it near a heating appliance or other high-temperature place.
- Do not allow cigarettes and other open flames to come close to the massage chair.


Massage chair installation



—Massage chair installation

How to move

- When transporting, please turn off the power switch and pull out the power plug from the socket.
- Push the back of the chair down first so that the rear wheels touch the ground, and then move it.
- To avoid damaging the floor, please place a mat or other material on the floor before moving the massage chair.

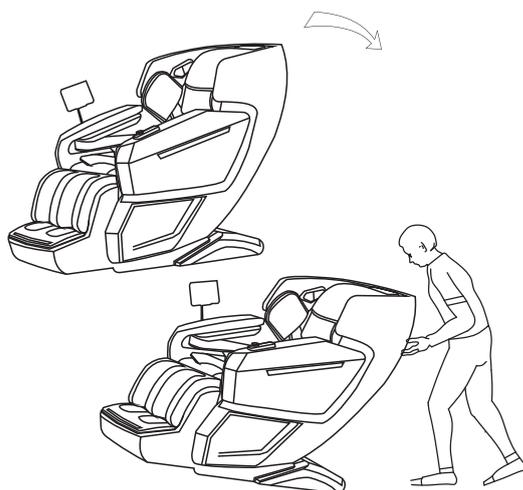

Notice

- If the aisle is narrower than 90cm, the massage chair cannot be transported normally.
- If you cannot transport the massage chair normally, please contact the dealer from whom you purchased it.
- Since the massage chair is heavy, be careful not to drop it when moving it in.
- Do not move the massage chair with people in it.
- Please unplug the power plug from the socket before moving the massage chair. (Otherwise, the power plug may be damaged, causing fire or electric shock.)




Installation of massage chair

Massage chair installation


Installation method

Massage chair unpacking

Lift the outer box up and remove the fixing post, take out the manual, power cords and other accessories.

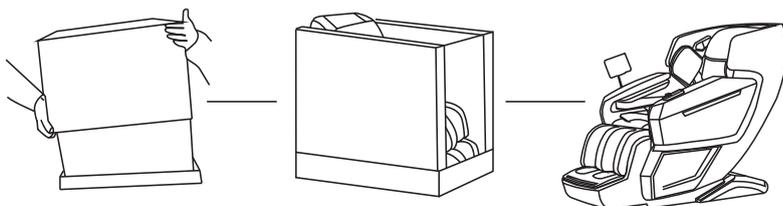

After unpacking the box, look for the main control switch under the back of the massage chair, connect the power cord, turn on the switch, and the massage is ready.

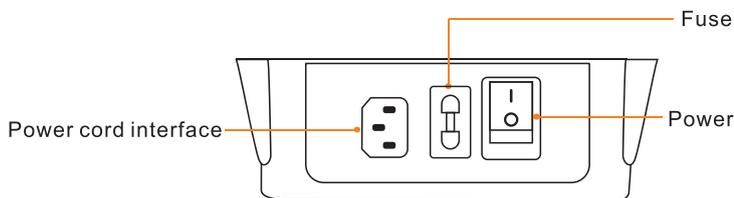

After the accessories are installed, insert the power plug into the 220V socket, turn on the power switch, and enable the hand controller to use it.

Maintenance

To clean dirt from leather and plastic parts, please follow the steps below

1. Turn off the power and pull out the power plug.
2. Soak a soft dry cloth in water or hot water containing 3% to 5% neutral detergent, and then wring it out thoroughly.
3. Wipe the dirt on the leather surface and plastic surface.
4. Use a dry cloth to wipe again.
5. Let it dry naturally.

Precautions

- Do not use alcohol, thinner, gasoline and other chemicals for cleaning
- Do not use a hair dryer to dry the surface.
- Never iron the fabric of the product.
- When performing maintenance, be sure to unplug the power plug and do not pull or plug it in with wet hands.
- If it will not be used for a long time, please use a dust cover (self prepare) to avoid being covered with dust.



Troubleshooting

Troubleshooting



Troubleshooting

If there is still an abnormality after performing the following inspections, to prevent accidents, please stop using it immediately, unplug the power plug from the socket, and contact the dealer from which you purchased it or the nearby customer service center.

Problem

Q1: The massage wheel stops during the use of the machine (the massage chair makes a "beep" sound, all the display icons and buttons on the remote control flash and then go out)

Q2:

Can't tilt

The calf cannot be raised and lowered (a "beep" sound is emitted, and the operation indicator of the remote control is turned off and stops running)

Q3: Massage does not reach the position of the shoulder or neck

Q4: This machine cannot run at all

- upper body (neck to waist) massage
- Calf and foot massage

Q5: Even if the warm function is turned off, the temperature will not come down.

Q6: The machine has been damaged

Q7: Power cord or power plug is abnormally hot

Cause and treatment

A1:

During use, when the massage wheel encounters excessive external force, the massage wheel may stop moving for safety reasons.

If all the markers and buttons on the controller begin to flash, this is extremely unusual.

Please turn off the power switch, wait about 10 seconds, then turn the power back on.

If the massage wheel stops again, adjust your sitting position, don't let your back completely on the backrest.

(Please pay attention to the relevant precautions, especially when using people with a weight of 100kg or more)

A2: If the unit encounters an obstacle or is too aggressive when using the unit, the unit will stop working for safety reasons.

When all the displayed contents and buttons are flashing, and the backrest or calf is stopped, turn off the power of the unit once, and then "on" the power again after about 10 seconds, then operate the unit.

A3: This situation is purely natural because it is equipped with a massage movement that can perform the massage action of the massage wheel. Not a malfunction of the machine.

A4: Power cord connection has been disconnected.

The power switch of this unit is not turned on.

No program selected after pressing the power button on the remote control

A5: Due to the structure of the heating bag, it will still feel warm due to the residual heat after heating.

Even if you turn off the heating function, when continue to massage the same part, it will make the body temperature rise.

A6-A7: In order to prevent accidents, please contact the after-sales center



Product parameters and specifications

Technical specifications

Product parameters and specifications

Product model	SATORI ELITE
Rated voltage	110V-240V~
Rated frequency	50Hz-60Hz
Rated power	150W
Working hours	The default is 15min, adjustable to 10-15-20-25-30min
Executive standard	GB4706. 1-2005 GB4706. 10-2008
Safety category	I class
Reset size	1460 × 760 × 1180 mm (Length * width * height)
Lie down for maximum size	1800 × 760 × 820mm (Length * width * height)

VOICE COMMANDS

	The command word	reply
1	HI-MASSAGE-CHAIR	HI
2	START-MASSAGE	OK
3	END-MASSAGE	OK
4	INCREASE-MASSAGE	OK
5	REDUCE-MASSAGE	OK
6	START-MODE-ONE	OK
7	START-MODE-TWO	OK
8	START-MODE-THREE	OK
9	START-MODE-FOUR	OK
10	START-MODE-FIVE	OK
11	START-MODE-SIX	OK
12	START-HEAT	OK
13	END-HEAT	OK
14	INCREASE-AIR	OK
15	REDUCE-AIR	OK
16	START-FOOT	OK
17	END-FOOT	OK
18	BACK-RISE	OK
19	BACK-DOWN	OK
20	BACK-STOP	OK
21	START-AIR	OK
22	END-AIR	OK
23	START-ZERO	OK
24	END-ZERO	OK